



FLUID METERING, INC.

Precision Valveless Metering Pumps and Dispensers

Over 59 Years of Precision Valveless Fluid Control & Over One Million Pumps




Medical, Analytical & Biotech Instrumentation
2018 Products Catalog

www.FluidMetering.com / 516 922 6050




One Dispenser / Pump for All OEM Applications

Valveless Syringing Aspirate & Dispense

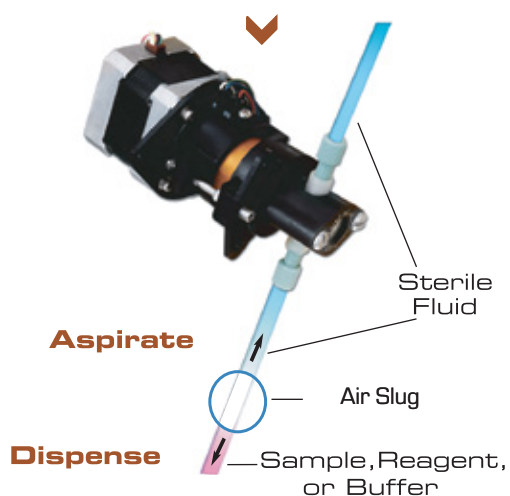

Fast Prime Flush & Wash

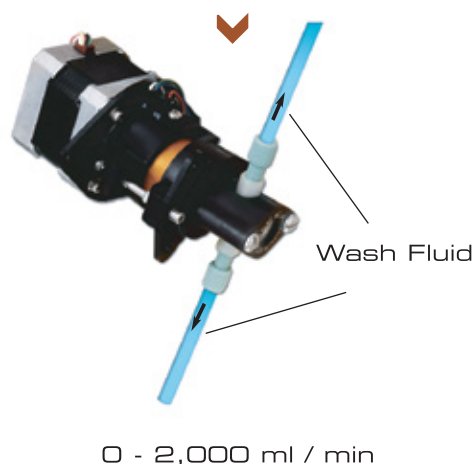


Continuous Dispensing



Continuous Metering



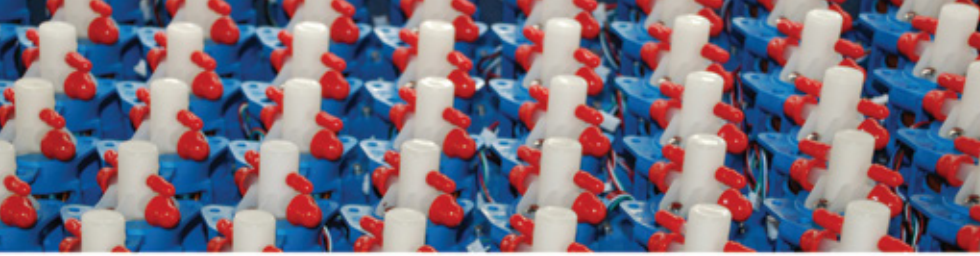


Why FMI?


Ultra-Precise Fluid Control. . . from Sub-Microliters to Liters

- Patented “No-Valve” Design**
 Eliminates problems and errors caused by valves which clog, leak, hang up, and require service.
- One Moving Part!**
 CeramPump® design utilizes a single, dimensionally stable, chemically inert CERAMIC piston and cylinder ensuring long term, drift free fluid control.
- Proven Performance!**
 Over 59 years OEM experience and more than a million pumps sold.
- Accuracy, Precision, & Reliability**
 Better than $\pm 1\%$ Measured in millions of “trouble-free” cycles.


Precision
Repeatability & degree of variation of a set of values

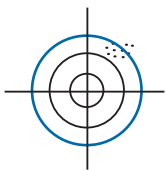

Accuracy
How close the average value is to the true value

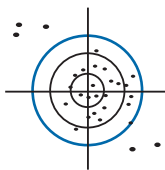

FMI Pumps
Accurate & Precise

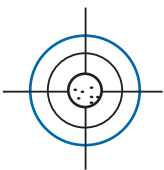




How Our Pump Works

The valveless pumping function is accomplished by the synchronous rotation and reciprocation of a ceramic piston in a precisely mated ceramic cylinder liner.

One complete piston revolution is required for each suction/discharge cycle.

Check out our web site at www.fluidmetering.com for an animation of our unique pumping principle.

Accuracy, Precision, & Reliability. . . FMI Pump Typical Flow Data

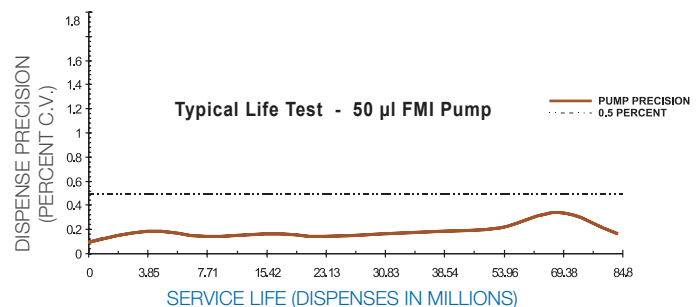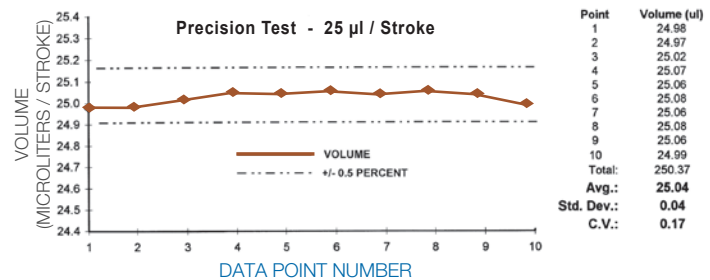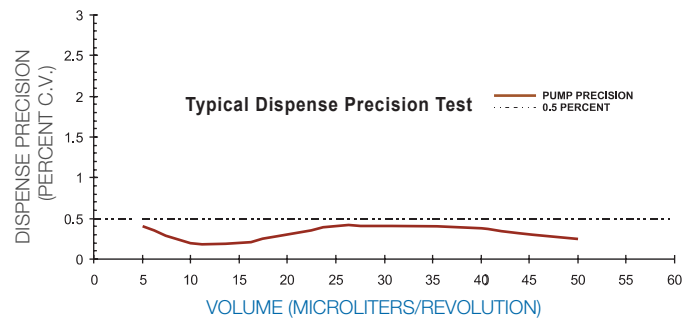

Celebrating 59 Years of Precision Valveless Fluid Control

Capabilities

- From miniature OEM Stepper Pumps to High Torque Hazardous Duty
- From 3 microliters per dispense up to 4.6 liters/min continuous flow
- Operating Temperatures Up to 350° F (158° C)
- Pressures up to 200 psig (13.8 bar)
- Over a Million OEM Pumps Sold

“STRH” Precision Adjustment Stepper Pump

Ideal for Prototyping & Applications
Requiring Frequent, Accurate Changes
in Dispense Volume

- Precision RH adjustable pump with integral stepper motor & opto sensor.
- Pump head displacement is variable using an easy-grip flow ring ring.
- Low dead-volume design ensures maximum bubble clearing.
- Sample, Dispense, Aspirate, Flush.



“STQP” Adjustable High Flow Stepper Pump

Ideal for Prototyping

- Precision, variable displacement “Q” Pump with integral stepper motor.
- Accommodates all “Q” style pump heads and RH pump heads (with RH/Q adaptor).
- Ideal for OEM applications where accurate & frequent displacement changes are expected.



“STH” Low Flow “STQ” High Flow OEM Pumps & Dispensers

Ideal for Medical, Analytical, & OEM Instrumentation Applications

- Isolation Gland now available for both STH & STQ models for handling crystallizing fluids.
- Duplex models available for dual channel or proportional metering & dispensing.
- Compact design.
- Factory calibrated.


“STH” Low Flow

Complete Pump Assembly	Max. Dispense Range Microliters / Revolution
STH00CKCLF	0 - 25 µl
STH0CKCLF	0 - 50 µl
STH1CKCLF	0 - 100 µl
STH2CKC	0 - 200 µl



“STQ” High Flow

Complete Pump Assembly	Max. Dispense Range Milliliters / Revolution
STQ1CKC	0 - .32 mL
STQ2CKC	0 - .72 mL
STQ3CKC	0 - 1.28 mL

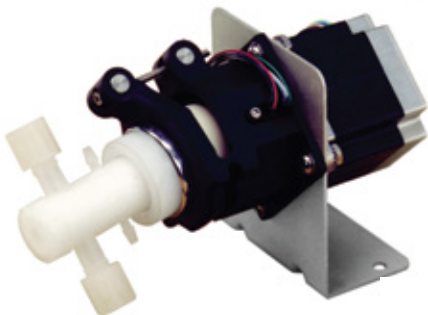

“RH” Miniature Pedestal Pump

Panel Mount high precision

- Miniature Motorless “RH” Pumps - Low Flow, High Precision.
- Ceramic & PVDF Fluorocarbon standard wetted materials.
- 0 - 100 Microliters per stroke.
- Precision 0.5% or better.
- Needs only 17 inch ounces of torque.
- Requires only 2 1/4" panel space.
- Easy-grip flow control ring.



“QP” Motorless Pedestal Pump

Pedestal Design Used Extensively in Industrial, Laboratory, & OEM Applications

- Easily handles slurries, suspensions, emulsions, solvents, viscous concentrates, gases, and more...
- Typically driven by belt, chain, or shaft coupling connected to your special motor drive.
- Minimal torque requirement of 35 inch ounces.



"STF" Fixed Displacement Pump

Ideal for OEM Instrumentation & Equipment for Medical, Analytical, & Industrial Applications

- Economical design with fixed displacement link.
- Displacement link can be customized for individual requirements.
- Precision stepper motors with opto sensors.
- Isolation gland available for crystallizing fluids.




"Multi Channel" Pump



"3-Way" Pump

Ideal for Low Volume Blending or Proportioning in Medical Diagnostic Instrumentation

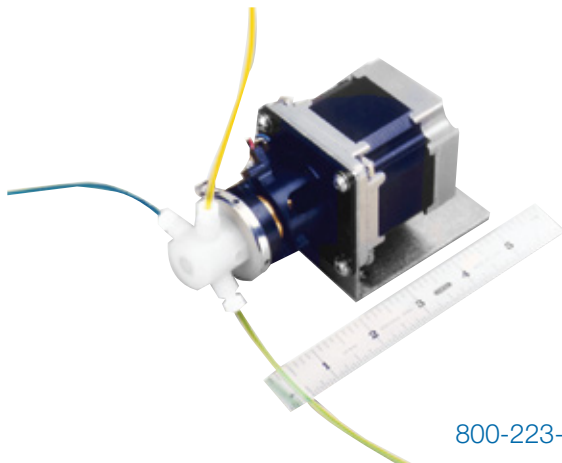

Brushless DC Pumps

Ideal for Waste/Transfer Applications in Medical, Analytical & Process Instrumentation

- Compact design with integral electronics.
- 24 V Brushless DC motor.
- Factory calibrated for flow rates from 20 mL/min to 400 mL/min.



DC Pumps

"The Standard" for Agricultural Spraying, Mosquito Control and Hydrogen Fuel Cell


- Micrometer-like flow adjustment to 260 mL/min.
- Extended motor shaft can be used for user supplied tachometer.



"ICST-02" Stepper Control

Programmable Control for All FMI Stepper Pumps

- Multiple Input and Output Connections.
- RS 232 Serial Port for PC Connection.
- MS Windows® Programming Software Included.
- Compact Size: 2.0" x 3.1" x 1.6" high (51 x 79 x 41 mm).

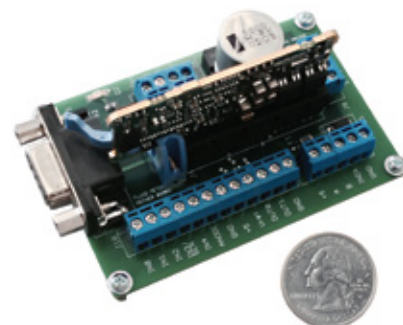

“Intelligent” Programmable Pump


Quick Start Control for All FMI Stepper Pumps

- FMI's STH Stepper Pump with integral programmable driver.
- Driver provides servo control of a stepper pump.
- 5 Programmable inputs, 2 programmable outputs.
- Multiple programming platforms including Visual Basic, C/C++, Delphi, Lab VIEW.
- Analog 0-5V, RS-232 serial, CANopen protocol supported.
- Resonance-free, quiet operation.
- EtherCAT (with optional module).



“400 μ L” Pump

- Compact Size - Ideal for OEM Instrumentation
- Nine (9) Stepper Drive Combinations provide a combined dispense range of 200 to 400 μ L and continuous metering up to 400 mL/min.
- Adjustable Link provides $\pm 25 \mu$ L displacement adjustment for each drive combination
- Four (4) Pump Head configurations available for maximum chemical compatibility, as well as to accommodate both low volume dispensing and high flow metering.

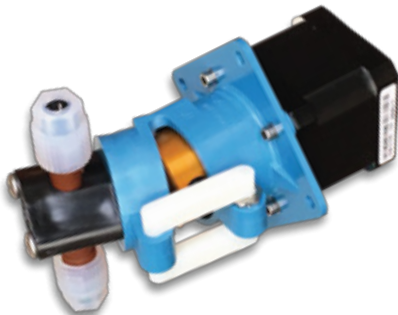

Duplex Ratio:Matic® Pumps


Ideal for Metering reagents and dilution of buffers in OEM medical, laboratory, and environmental analyzers. Production fluid dispensing for manufacture of contact lenses, button cell batteries, and reagent test kits.

- Two valveless pumps controlled by a single variable speed drive.
- Pump displacements independently variable for proportional metering.
- Vary flow rate while maintaining ratios of fluid dispensed or metered.
- Stepper motor and variable DC drive configurations.
- Up to 500:1 Ratio.


“STH/Q”



“ST2RH”

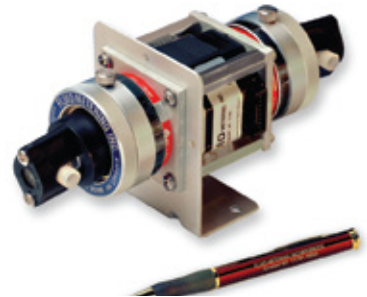

“ST2Q”

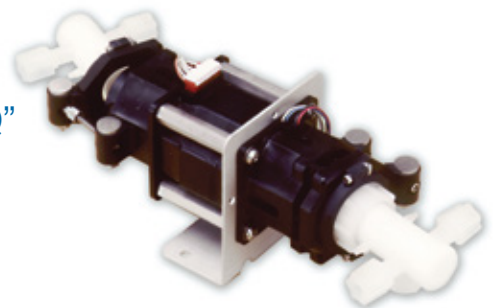

“ST2H”



Who We Are

FMI pioneered the first patented valveless rotating and reciprocating piston metering pump concept and has been delivering pumping excellence and precise fluid control for over 59 years.

Quality

We take quality seriously and back it up, not only as an ISO 9001:2008 facility, but far beyond! Most products are WEEE & RoHS Compliant.

Our Mission Statement . . .
100% Quality, 100% On-Time Delivery

. . . is supported by our valued
OEM supplier awards.



Engineering Design & Development

Our Engineering Team incorporates over 59 years of OEM design experience to meet specific customer & application requirements. With the knowledge and the necessary tools, our engineers have developed the most precise and reliable valveless dispensing and metering pumps available.

eSupport (FMI web site)

Need product and technical information immediately? Check our web site at www.fluidmetering.com and have instant access to product specifications, application information, literature downloads, and an animation of our unique CeramPump® valveless pumping principle.

Also featured in our web site is LiveHelp, which provides a one on one connection between our customers and FMI's application specialists.



Have questions?
Chat live with an FMI
application specialist at
www.fluidmetering.com

Applications. . . We Do It All !

MEDICAL

For precise dispensing, aspirating, rinsing, mixing systems and syringe pump replacement in diagnostic, clinical chemistry and medical equipment manufacturing.

SPRAYING SYSTEMS

For injection of insecticides, herbicides, and agricultural nutrients, as well as for ULV spray equipment.

ENVIRONMENTAL & POLLUTION CONTROL

For sampling stack gases, ground water & wastewater, as well as injection of monomers, polymers, and chemicals for water & waste treatment, TCLP and more. . .

OFFICIAL UK DISTRIBUTOR:

Michael Smith Engineers Limited

Website: www.michael-smith-engineers.co.uk

Freephone: 0800 316 7891



INSTRUMENTATION

For all kinds of precision instruments and monitors including titrators, TOC, S₂O monitors, chromatographic systems & humidity control.

DISPENSING SYSTEMS

For dispensing of solvents, UV adhesives, lubricants, reagents, and mercury in the manufacture of electronics, pharmaceuticals, medical disposables, computer hard drives, audio speakers, and calibration equipment.

FOOD & DAIRY

For candy coating and polishing, as well as, vitamin fortification, addition of flavors, colors, preservatives and a variety of other ingredients used to enhance food products.



5 Aerial Way, Suite 500
Syosset, NY 11791 USA
Toll free: 800-223-3388
Tel: 516- 922-6050
Fax: 516-624-8261
Email: pumps@fmipump.com